

Presenting the
SPECTROMAX^x LMX10
ARC/SPARK OES ANALYZER



**Analytical improvements and more
for next-generation material control**

Best-selling analysis to optimize material control

SPECTROMAXx

Metal producing and fabricating plants require truly advanced elemental analysis at every step — from incoming materials to in-process testing to final inspection for outgoing quality. One analytical instrument measures up every time.

With outstanding repeatability, reproducibility, and reliability, SPECTROMAXx is perhaps the industry's best-selling arc/spark optical emission spectrometry (OES) analyzer. Its fast, accurate, cost-effective measurements add certainty to critical supply chains.

Users get ultrafast information on changing process conditions. Also drastically reduced cost of ownership — with lower consumables, plus advanced diagnostics and easy maintenance to prevent expensive downtime. And this tenth-generation SPECTROMAXx responds to evolving plant requirements with a host of new improvements, as shown below.



FAST, SIMPLE STANDARDIZATION WITH ICAL 2.0

Using conventional analyzers, standardization can take 30+ minutes, demand multiple samples, and require reruns whenever site conditions change. By contrast, in most cases SPECTRO's proprietary iCAL 2.0 calibration logic needs only 5 minutes and a single sample per day. Plus it automatically compensates for most changes in environmental temperature or pressure. So SPECTROMAXx delivers outstanding stability, productivity, and savings, time after time.

NEW ANALYTICAL PERFORMANCE

SPECTROMAXx LMX10 features completely new factory calibrations, with extended calibration ranges and element selections; optimized source excitation parameters; and improved limits of detection (LODs). Example: for iron (Fe), aluminum (Al), and copper (Cu) matrices, the analyzer surpasses its already impressive performance with 30% to 40% better LODs.

NEW ARGON COST SAVINGS

Previous SPECTROMAXx models significantly decreased their consumption of expensive argon (Ar) gas. A brand-new feature now allows total shutdown of argon flow during standby, for periods configurable by the operator. Users save more consumables costs than ever!

NEW EASE OF USE

SPECTROMAXx now provides routine, trouble-free analysis of 10 matrices, 65 methods, and 59 elements — via convenient controls for operation; easy access for use and maintenance; and enhanced software features such as quick-check programs, virtual type standards, and spectrum scans.

SPECTROMAXx

reliability for incoming and outgoing material control

SPECTROMAXx AND SPECTROMAXx ADVANCED


The new generation offers two models, differing only in their optical systems.

Basic **SPECTROMAXx** features a single air optic, with state-of-the-art CMOS sensors proven in the flagship SPECTROLAB S. It extends the relevant and applicable elemental wavelength range from 233 nanometers (nm) to 670 nm. A new option analyzes relevant wavelengths to 766 nm, covering potassium (K).

SPECTROMAXx Advanced utilizes the same high-resolution CMOS technology, with the air optic covering the same wavelengths as the basic model, from 233 nm to 670 nm. But an added UV optic handles a lower spectral range, from 120 nm to 235 nm. This includes a new capability to analyze oxygen (O) in copper. And again, a new option analyzes relevant wavelengths to 766 nm, covering potassium.



HIGH-PERFORMANCE HARDWARE

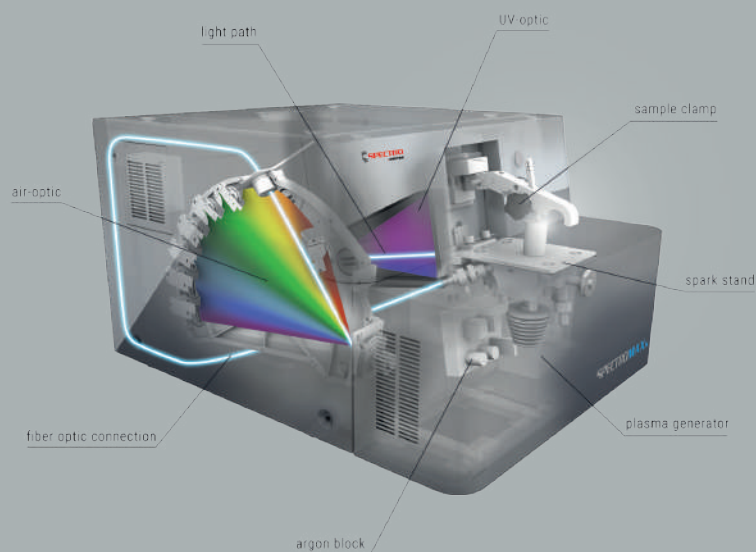

HARDWARE NOTE:

The latest SPECTROMAXx eliminates a magnesium fluoride (MgF_2) window at the spark stand — saving costs and maintenance replacement efforts.

EASY-TO-USE SOFTWARE

Even for less experienced personnel, SPECTRO's SPARK ANALYZER Pro software provides effortless operation. Favorites include application profiles tailored to preset user requirements; automatic program selection of sub-methods for given materials; spectra storage for later recalculations; multiple data export formats; and software-only upgrades for new lines, methods, or matrices.

Newest software highlights:

- **Backup tool** — for easy backups to network, other drives, etc.
- **Virtual type standards** — to improve accuracy by combining optimal element samples for type correction
- **Extended grade library** — for iron, aluminum, copper
- **Quick-check programs** — to identify & analyze main elements in under 12 seconds: for iron & aluminum materials, & now including analysis of copper
- **Spectrum scan (optional)** — to compare spectral lines of multiple samples in the same view
- **Shift change protocol (optional)** — to summarize all activities in a given shift: for handover to next operator, or QA optimization

SPECTROMAXx additional advantages



ADAPTERS FOR VARYING NEEDS

Material control presents a wide range of analysis requirements. Two specific challenges: adjusting to samples' differing shapes and sizes, and optimizing positioning on the spark stand. Fortunately, SPECTROMAXx's adapter kit offers a variety of flexible, easy-to-use solutions.



Small parts adapter kit



Adapter for wires with a diameter of 0.8-3 mm



Adapter for wires/screws or small parts with a diameter of 0.8-10 mm



Adapter for small pieces



SUPERIOR SERVICE AND SUPPORT

AMECARE Performance Services maximize uptime for all the world-class elemental analyzer products and services from SPECTRO Analytical and associated companies. The program is staffed by hundreds of experienced service engineers in 50 countries. They provide high-value, customized support designed to ensure the longest possible equipment life— plus optimum performance for high-volume, high-productivity analysis.

Ask about AMECARE virtual or on-site demos, SPECTRO PROTEKT secure global remote monitoring, proactive performance maintenance, performance upgrades, applications solutions, consultation, targeted training, and ongoing support.

COMPLETE LINE OF METAL ANALYZERS

This next-generation SPECTROMAXx model takes its place among today's most comprehensive suite of advanced stationary metal analyzers. It stands beside SPECTRO's flagship SPECTROLAB S, a market leader for ultra-precise analysis, plus the entry-level SPECTROCHECK, offering high performance and dependability at a low, very competitive price.

Complementing these are SPECTRO's mobile metal analyzers: the SPECTROTEST mobile analyzer, the SPECTROPORT portable metal analyzer, and the SPECTRO xSORT XRF handheld analyzer.

www.spectro.com

GERMANY

SPECTRO Analytical Instruments GmbH
Boschstrasse 10
D-47533 Kleve
Tel. +49.2821.892.0
spectro.sales@ametek.com



U.S.A.

SPECTRO Analytical Instruments Inc.
50 Fordham Rd
Wilmington 01887, MA
Tel. +1 800 548 5809
+1 201 642 3000
spectro-usa.sales@ametek.com

QS Instruments Sdn. Bhd.

No. 2, Jalan Jururancang U1/21,
Hicom Glenmarie Industrial Park,
Seksyen U1, 40150 Shah Alam,
Selangor, West Malaysia
Tel: +603-5882 6866
Fax : +603-5567 0813
spi-sales@qesnet.com
www.qesnet.com

Subsidiaries:

- ▶ **FRANCE:** Tel. +33.1.3068.8970, spectro-france.sales@ametek.com
- ▶ **GREAT BRITAIN:** Tel. +44.1162.462.950, spectro-uk.sales@ametek.com
- ▶ **INDIA:** Tel. +91.22.6196.8200, sales.spectroindia@ametek.com
- ▶ **ITALY:** Tel. +39.02.94693.1, spectro-italy.sales@ametek.com
- ▶ **JAPAN:** Tel. +81.3.6809.2405, spectro-japan.info@ametek.co.jp
- ▶ **SOUTH AFRICA:** Tel. +27.11.979.4241, spectro-za.sales@ametek.com

SPECTRO operates worldwide and is present in more than 50 countries. For SPECTRO near you, please visit www.spectro.com/worldwide

© 2022 AMETEK Inc., all rights reserved, subject to technical modifications • F-22, Rev. 0 • 80902232 • Photos: SPECTRO, Adobe Stock
Registered trademarks of SPECTRO Analytical Instruments GmbH • **SPECTRO**: USA (3,645,267); EU (005673694); "SPECTRO": EU (009693763);
„SPECTROTEST": EU (004206173), USA (4,103,718); "SPECTROPORT": EU (1283844), USA (5,056,898), China (1283844); xSORT: EU (7058456), USA (3,767,555);
„SPECTROLAB": EU (1069339), USA (4,103,747)